



PROACTIV®



NJ1 adaptive bike

THE NJ1 ADAPTIVE BIKE – FLEXIBILITY MEETS COMFORT

The **NJ1 adaptive bike** is a lightweight and sturdy manual propulsion aid for leisure and sport. The bike simply attaches to the wheelchair, enabling ergonomic movement and allowing users to cover longer distances. The larger drive wheel makes it easier to negotiate uneven surfaces and obstacles, whilst high-quality gear and brake components ensure precise control. Thanks to its versatile features, the NJ1 can be individually configured to suit everything from touring to sport cycling.



#LightweightAndPowerful

The NJ1 adaptive bike combines low weight with high stability and excellent handling. This makes it the ideal companion for leisure, sport and longer rides.



#EasyHandling

The handy parking stand makes it easier to attach and detach the trailer and ensures it is parked safely. It can be folded away to save space whilst driving.



#StableAdaptation

With a central or front adapter, the NJ1 adaptive bike can be quickly and securely attached to a wheelchair – for reliable adaptation and easy handling.



#YourNJ1YourChoice

Many options for your perfect handbike. Every NJ1 bike is customised. A wide range of options and accessories make the NJ1 your bike.



WIDE RANGE OF CONFIGURATION OPTIONS FOR MAXIMUM CUSTOMISATION.

The **NJ1 adaptive bike** offers flexible configurations for an optimised riding experience. Extensive equipment options, such as a practical kickstand, and the choice of various high-quality gears and brakes make handling comfortable and easy.

SELECTABLE OPTIONS FROM PRO ACTIV:




Individual handlebar widths and crank lengths, as well as the adjustable bottom bracket support, ensure the bike is perfectly tailored to your needs.



Perfect gearing. The NJ1 adaptive bike offers a wide selection of gears, from derailleurs gears to bottom bracket gears. Whether mechanical or electric – you decide!



Reliable brakes, maximum control. With disc brakes, rim brakes or the optional back-pedal brake, you're in safe hands.

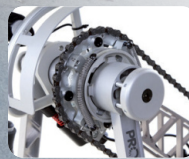

The optional back-pedal brake with crank release function allows the cranks to be disengaged, enabling the wheelchair to be moved backwards using the handrims – ideal for manoeuvring in confined spaces.

TECHNICAL SPECIFICATIONS

Drive type:	manual with gear shift
Special features:	stable with excellent handling, with a range of configurations from touring to sports bikes and off-road
Total weight without adapter:	starting at 8.3 kg
Pack size:	starting at 115 x 43 x 50.8 cm
Wheel sizes:	20" 24"
Max. standard load weight:	up to 140 kg trailer load and 10 kg load capacity
Max. load weight for Heavy Duty:	up to 170 kg towing capacity; higher towing capacities available on request
Gear systems:	derailleur gears or hub gears with operation on the handle bottom bracket gears
Brakes:	disc brake rim brake back-pedal brake

FOLLOW US ON SOCIAL MEDIA



PRO ACTIV

PRO ACTIV Reha-Technik GmbH

Im Hofstätt 11
D-72359 Dotternhausen

info@proactiv-gmbh.de

www.proactiv-gmbh.com