



PROACTIV®



SPEEDY

F3

MAXIMUM PERFORMANCE FOR EXPERIENCED USERS.

The SPEEDY F3 is a highly stable rigid-frame wheelchair for experienced users. Its fully welded frame ensures outstanding handling on both flat and uneven terrain. The seating position and tilt point are individually adjusted to ensure maximum comfort and optimum power transfer. Thanks to 14 selectable frame parameters, the SPEEDY F3 can be perfectly tailored to individual needs.



#Reliability

The SPEEDY F3 impresses with a sturdy and durable frame that ensures maximum mobility and is compatible with various auxiliary drive systems and adaptations.



#MaximumPerformance

The SPEEDY F3 offers unbeatable handling – on both flat and uneven terrain. The seating position and tilt point are individually adjusted to suit the user. Ideal for experienced wheelchair users.



##LightweightAndStable

The fully welded frame offers maximum stability whilst keeping weight to a minimum. A customised frame extension with frame retraction and an adjusted knee angle enables precise leg movement.



#Independence

Thanks to the easily foldable backrest and toolfree removal of assemblies, the SPEEDY F3 is easy to transport.



MAXIMUM FLEXIBILITY.

Your **SPEEDY F3** is manufactured precisely to your requirements. You can choose from a wide range of options – from frame sizes, footrests, brakes and push handles to backrest and seating systems, clothing guards, armrests and drive and steering wheels. Your wheelchair, your choice.

UNIQUE, CUSTOMISABLE OPTIONS FROM PRO ACTIV:



The 7-way angle-adjustable and foldable backrest offers maximum flexibility for the ideal sitting position.



Three seating systems for maximum comfort: **Open strap system** for individually adjustable pressure distribution. **Body Contour** with cushioning effect and pressure ulcer prevention. **Seat panel** for stable seating without sagging.



You can choose between a permanently welded or a removable wheelbase extension for greater traction when using a handbike or towing device.






With the Heavy Duty equipment, your **SPEEDY F3** is designed for extreme loads and is suitable for user weights up to 200 kg.

TECHNICAL SPECIFICATIONS

Seat width:	25 - 55 cm (in 1 cm steps)
Seat depth:	23 - 55 cm (in 1 cm steps)
Back height:	20 - 60 cm (in 1 cm steps)
Lower leg length:	35 - 52 cm (in 1 cm steps)
Rear seat height:	35 - 50 cm (in 1 cm steps)
Front seat height:	36 - 58 cm (in 1 cm steps)
Tipping point:	-5 - 17 cm (in 1 cm steps)
Total product weight:	Starting at 8.5 kg
Max. standard load weight:	Up to 120 kg
Max. load weight for Heavy Duty:	Up to 200 kg appropriate individual dimensions upon request

Data subject to design changes.

FOLLOW US ON SOCIAL MEDIA

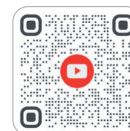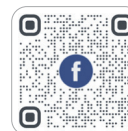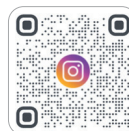

PRO ACTIV®

PRO ACTIV Reha-Technik GmbH

Im Hofstätt 11
D-72359 Dotternhausen
info@proactiv-gmbh.de

www.proactiv-gmbh.com

Rev. 07/2026
Subject to technical modifications, printing errors & colour deviations.