



PRO[®]ACTIV



HUSK-E


POWERFUL AND INCREDIBLY SMALL DURING TRANSPORT

The **HUSK-E** adaptive bike combines flexibility, performance and comfort in a compact design. The foldable bottom bracket support ensures a small pack size – ideal for transport. There are three versions to choose from: the powerful Shimano model for mountainous terrain, the quiet and energy-saving Neodrives version for flatter terrain and the manual version for sporty self-propulsion. For older children and teenagers, the HUSK-E is available with a short bottom bracket support. A customised gear selection, ergonomic adjustment options and efficient brakes ensure optimum riding comfort and safety.



#PerformanceThatInspires

Powerful assistance up to 25 km/h. Whether flat or hilly: The drive system of the HUSK-E supports you up to 25 km/h. Cover longer distances effortlessly and enjoy the freedom of movement.

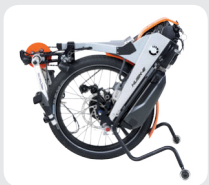

#CompactAndFlexible

The HUSK-E revolutionises the transport of handbikes: the fold-down bottom bracket support allows the adaptive bike to be quickly and easily reduced in size. Ideal for travelling and space-saving storage.



#EnduranceForEveryTour

Large battery capacity and dual battery option. With ranges of up to 85 km and the option of a double battery holder, even long tours are no problem.



#ThreeVariantsForEveryRequirement

HUSK-E is available as a Shimano (powerful for mountainous terrain), Neodrives (quiet, energy-saving for flatter terrain) or as a manual version – depending on your needs.



FLEXIBLY CONFIGURABLE FOR EVERY NEED.

The **HUSK-E** adapts to you – whether with a Shimano or Neodrives drive system or as a manual bike without a motor. Various gears, brakes and ergonomic adjustments ensure maximum comfort. The fold-down bottom bracket support, practical kickstand and long reach make it suitable for touring and easy to transport.

UNIQUE, SELECTABLE OPTIONS FROM PRO ACTIV:




Reliable adaptation and easy handling. The HUSK-E is quickly attached to the wheelchair using the centre or front adapter. The practical parking stand also makes coupling and uncoupling easier.



Perfect ergonomics for maximum comfort. Individual grip widths and crank lengths as well as the adjustable bottom bracket support make the HUSK-E your ideal companion.



The manual version of the HUSK-E offers a sporty riding experience – without a motor. Ideal for anyone who wants to ride actively and improve their endurance.

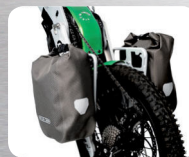



Luggage bags including holders are available as an option. The bags with a volume of 5.8 litres offer enough space for your shopping.

TECHNICAL SPECIFICATIONS

Motors: [if equipped with electric support]	Mid-mounted motor: Shimano 85 Nm 250 W 2.7 kg
	Hub motor: Neodrives 12/40 Nm 250 W 4.6 kg
Batteries: [if equipped with electric support]	Shimano 36 V 14 Ah 504 Wh 2.6 kg
	Range approx. 45 - 75 km
	Shimano 36 V 17.5 Ah 630 Wh 3.5 kg
	Range approx. 55 - 85 km
	Neodrives 36 V 16 Ah 576 Wh 3.6 kg
	Range approx. 45 - 75 km
Total weight without adapter:	For electric version: from 14.1 kg (without battery) For manual version: from 11.4 kg
Pack size:	Starting at 80 x 43 x 80 cm
Wheel sizes:	20"
Pushing aid/Starting aid: [if equipped with electric support]	Up to 4 km/h forwards (Shimano) Up to 6 km/h forwards and up to 1.5 km/h in reverse (Neodrives)

Data subject to design changes.

FOLLOW US ON SOCIAL MEDIA

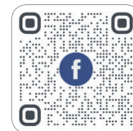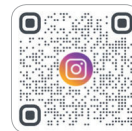

PROACTIV®

PRO ACTIV Reha-Technik GmbH

Im Hofstätt 11
72359 Dotternhausen
info@proactiv-gmbh.de

www.proactiv-gmbh.com

Rev. 05/2025

Subject to technical modifications, printing errors & colour deviations.