



PRO^{ACTIV}



30cm



LIFT activ


WITH THE LIFT ACTIV YOU CAN REACH MORE.

The **LIFT activ** with seat height adjustment makes you flexible in everyday life as well as in different professions and workplaces. You can effortlessly reach machines, furniture, kitchen fittings, printers and shelves – often without the need for any modifications. Weighing just 14.5 kg, the LIFT activ impresses with its lightweight, dynamic design and ease of transport. As well as the mechanical version, the LIFT activ is also available in an electric version. This allows you to adjust the seat height comfortably and continuously at the touch of a button – ideal for even greater independence, efficiency and ergonomics in everyday life and at work.



#Independence

With the LIFT activ, you can master everyday life with even more self-determination. When shopping, you can reach goods from the upper shelves without outside help and are independent.

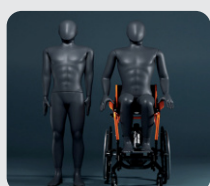

#AtEyeLevel

The infinitely variable seat height adjustment raises the seat by up to 30 cm. In this position, the wheelchair user is at eye level with a person standing about 1.75 m tall.




#ComfortAtTheTouchOfAButton

The LIFT activ electric allows you to adjust the seat height continuously at the touch of a button – precisely and effortlessly, for greater independence in everyday life.



#Independence

Thanks to the weight optimization, the foldable backrest and the tool-free removal of assemblies, you can easily transport your LIFT activ.

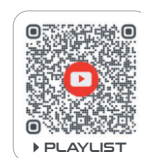

YOUR KEY TO MORE INDEPENDENCE.

Your **LIFT activ** will be manufactured individually for you, to suit your needs. You have a large selection of equipment – from footrests and brakes to push handles, seat and back systems as well as drive and caster wheels. So your LIFT activ is exactly as you need it.

UNIQUE, SELECTABLE OPTIONS FROM PRO ACTIV:



The 7-way angle-adjustable and foldable backrest offers maximum flexibility for the ideal sitting position.



The backrest, which can be folded down at the lumbar joint, allows an optimal pack size for transport – quickly and easily.



The removable and swivelling frame stem with split footrest offers plenty of space for your transfer.




The footrest, which can be folded towards the back, is easy to handle and flexible in everyday life.

TECHNICAL SPECIFICATIONS

Seat width:	35 - 46 cm (in 1 cm steps)
Seat depth:	35 - 48 cm (in 1 cm steps)
Back height:	20 - 48 cm (in 1 cm steps)
Lower leg length:	30 - 48 cm (in 1 cm steps)
Rear seat height:	41 - 43 cm (in 1 cm steps)
Front seat height:	45 - 51 cm (in 1 cm steps)
Tipping point:	6 - 10 cm (in 1 cm steps)
Total product weight:	mechanical: from 14.5 kg electric: from 23 kg
Max. standard load weight:	up to 100 kg

FOLLOW US ON SOCIAL MEDIA

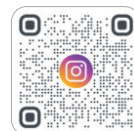

PRO ACTIV

PRO ACTIV Reha-Technik GmbH

Im Hofstätt 11
D-72359 Dotternhausen
info@proactiv-gmbh.de

www.proactiv-gmbh.com