



PROACTIV®



THE CHILDREN'S WHEELCHAIR FOR THE ADVENTURE OF LIFE.

The **BUDDY** is a compact, highly adaptable wheelchair for children and young people that grows with them and is ready for anything: everyday life, therapy, the playground or big outings. It promotes independence, mobility and therapeutic development, and supports children safely and reliably for many years to come. With its sturdy frame construction, intelligent growth-adjustment feature and numerous options, the BUDDY adapts perfectly to the child's life – not the other way round.



#GrowAlongFeature

The BUDDY grows with you – in width, height and depth. Thanks to its finely adjustable seat height and tilt-point adjustment, it adapts to your child's needs. This means the BUDDY remains a faithful companion for years to come.



#HealthyPosture

The BUDDY promotes a healthy sitting posture and encourages coordinated movements – the abducted frame design ensures additional legroom and optimal positioning options.



#CleverBraking

The BUDDY offers a range of knee-lever braking systems tailored to the needs of both the child and the accompanying adult. The integrated knee-lever braking system in the clothing guard is particularly child-friendly.



#Support

Two backrest systems ensure optimum seating comfort: a lightweight, adjustable backrest with individual adjustment options, or a sturdy backrest shell with lumbar support.




CUSTOMISABLE TO SUIT EVERY NEED

The **BUDDY** can be tailored precisely to your own requirements with a wide range of options. Whether it's a one-hand drive, therapy table, outdoor bike attachment, flexible footrests or various braking and back support systems – every accessory adapts to the child's individual needs.

UNIQUE, CUSTOMISABLE OPTIONS FROM PRO ACTIV:



The 'footrest that folds up towards the back' makes it easier to transfer in and out of the wheelchair for children who are still able to stand.

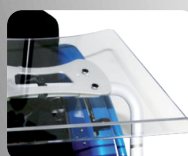

The therapy table offers a practical work surface, can be removed without tools and can be infinitely adjusted in height and angle for optimal adjustment.



With the **outdoor wheel attachment**, the BUDDY is suitable for off-road use on pavements, gravel and forest paths. The caster wheels are lifted off the ground and the 12" wheel overcomes bumps effortlessly. Thanks to the quick-release clamp, the wheel attachment can be attached easily and without tools.






TECHNICAL SPECIFICATIONS

Seat width:	20 - 46 cm (in 1 cm steps)
Seat depth:	22 - 40 cm (in 1 cm steps)
Back height:	20 - 48 cm (in 1 cm steps)
Lower leg length:	20 - 47 cm (in 1 cm steps)
Rear seat height:	31 - 37 cm [20" tyres] 33.5 - 39.5 cm [22" tyres] 36 - 42 cm [24" tyres] (in 1 cm steps)
Front seat height:	40 - 47 cm (in 1 cm steps)
Tipping point:	5 - 10 cm (in 1 cm steps)
Total product weight:	Starting at 9.5 kg
Max. standard load weight:	Up to 100 kg

Data subject to design changes.

FOLLOW US ON SOCIAL MEDIA

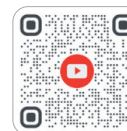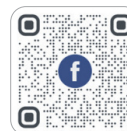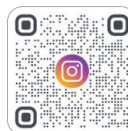

PROACTIV®

PRO ACTIV Reha-Technik GmbH

Im Hofstätt 11
D-72359 Dotternhausen
info@proactiv-gmbh.de

www.proactiv-gmbh.com

Rev. 06/2026

Subject to technical modifications, printing errors & colour deviations.