




PROACTIV®

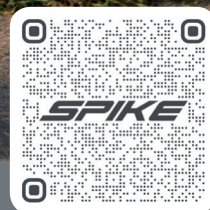

**SPIKE**

## COMPACT & LIGHTWEIGHT FOR MORE MOBILITY

The **SPIKE adaptive bike** is the ideal solution for children, teenagers and adults who are looking for a lightweight and compact bike. Thanks to its ergonomic movement sequences, it allows you to enjoy longer journeys without having to do without your own wheelchair. With its larger drive wheel, it masters obstacles and uneven surfaces with ease, while the braking system ensures safety on slopes. With a weight from 7.3 kg, a high-quality hub gear and a small turning radius, it is perfect for narrow alleyways, the school or work run and city tours.



### #CompactManoeuvrableSimple

The SPIKE adaptive bike impresses with its simple features and compact dimensions. Perfect for children, teenagers and adults who are looking for a lightweight bike for city tours, commuting to work or school.




### #LightweightAndTransportable

With a weight from 7.3 kg and a pack size of just 98 x 43 x 50.8 cm, the SPIKE is particularly easy to transport and fits easily into the car.



### #FlexibleCustomisable

Thanks to the versatile adjustment options on the bottom bracket support and the removable frame adapter, the SPIKE grows with your requirements.

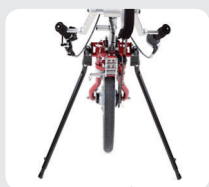

### #EasyHandling

The practical parking stand makes coupling and uncoupling easier and ensures safe parking. It can be folded up to save space while travelling.





## PERFECTLY ADAPTED – FOR EVERY RIDE.

The **SPIKE adaptive bike** offers you versatile configuration options for maximum comfort and safety. The grow-with-you function, various gears and effective brakes ensure your perfect riding experience. Thanks to its compact design and easy handling, it is ideal for everyday use and leisure.

### SELECTABLE OPTIONS FROM PRO ACTIV:



Reliable adaptation and easy handling. The SPIKE is quickly attached to the wheelchair using the centre or front adapter. The practical parking stand also makes coupling and uncoupling easier.



Perfect ergonomics for maximum comfort. Individual grip widths and crank lengths as well as an adjustable bottom bracket support make the SPIKE your ideal companion.



Effective brakes. The integrated coaster brake in the hub gear ensures simple and intuitive operation. In addition, a high-quality rim brake provides reliable braking power.





The numerous configuration options and extensive range of accessories allow you to personalise your equipment – for pure riding pleasure.

## TECHNICAL SPECIFICATIONS

Drive type:	Manual with gear shift or single speed
Special features:	Light   compact   simple   mostly used as a bike for children & young people or light adult bike
Total weight without adapter:	Starting at 7.3 kg
Pack size:	Starting at 98 x 43 x 50.8 cm
Wheel sizes:	16"   20"
Max. standard load weight:	Up to 75 kg trailer load and 5 kg load capacity
Gear systems:	Derailleur gears or hub gears with operation on the handle   bottom bracket gears
Brakes:	Disk brake   rim brake   back-pedal brake (with hub gears)

### FOLLOW US ON SOCIAL MEDIA

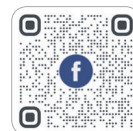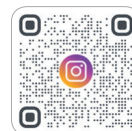

**PRO ACTIV®**

PRO ACTIV Reha-Technik GmbH

Im Hofstätt 11  
72359 Dotternhausen  
info@proactiv-gmbh.de

[www.proactiv-gmbh.com](http://www.proactiv-gmbh.com)