



PROACTIV®



WHEEL-E

MAXIMUM DRIVE POWER AND CONTROL FOR ANY TERRAIN

The **WHEEL-E** towing device transforms your wheelchair into a powerful companion for any route. With its 20" wheel, wide tyres and low centre of gravity, you can master paved paths, forest trails and inclines with ease. The battery offers a range of up to 45 km, extendable to 90 km with a double battery. Thanks to cruise control and stable steering geometry, you can enjoy maximum comfort and safety. Versatile accessories enable customisation – for independent mobility on any terrain!



#IndependentOnTheRoad

Activate the throttle with a flick of the wrist and roll effortlessly over pavements, forest paths or inclines thanks to the 20" wheel and wide tyres.



#MoreSpeed

Depending on the equipment, you can reach speeds of 6 km/h, 15 km/h (optionally with individual operating licence) or 24 km/h. Whether in the city, off-road or on private roads – the choice is yours.



#SteepIncline

Up to 15 % climbing ability. With the powerful Triride hub drive, you can even tackle challenging inclines. Perfectly combined with 5 support levels and a reverse gear for more flexibility.



#DoubleSoFar

Whether you're travelling short or long distances: up to 45 km range with one battery. With the optional double battery, you can enjoy double the range – and the interchangeable plug connection makes it easy to change the battery.



CUSTOMISED FOR YOU.

The **WHEEL-E** adapts perfectly to your needs: With an adjustable steering column in length, angle and height, you can find the ideal handlebar position. The flexible adapter can also be adjusted in height, length and angle – for maximum comfort and control.

UNIQUE, SELECTABLE OPTIONS FROM PRO ACTIV:



The sturdy centre or front adapter can be fitted in a few simple steps and hardly adds any weight to the wheelchair.



Cruise control for relaxed journeys. Enjoy constant speed at the touch of a button. No adjustment, no stress – just pure driving pleasure.



Individuality counts. With accessories such as mudguards, pannier racks and rear-view mirrors, the WHEEL-E can be customised to your needs – even for functionally impaired users.



150 kg are already possible as standard. A towing capacity of up to 170 kg is possible subject to a technical inspection and the fitting of a heavy-duty adapter. Higher towing capacities are available on request.






TECHNICAL SPECIFICATIONS

Motor:	Hub motor: Triride I max. 1000 W
Battery:	Triride 48 V 10.5 Ah 504 Wh Range approx. 25 - 45 km 3.3 kg
Double battery:	Possible; switch via change-over device on the bottom bracket support
Total weight without adapter:	Starting at 14.7 kg (without battery)
Max. standard load weight:	Up to 150 kg trailer load and 10 kg payload
Pack size:	Starting at 90-110 x 63 x 62 cm
Wheel sizes:	20"
Drive system speed:	6 km/h (Switzerland: 10 km/h), 15 km/h (only permitted with individual operating licence in Germany in public areas and in road traffic), 24 km/h (only permitted on private property in Germany)

Data subject to design changes.

FOLLOW US ON SOCIAL MEDIA

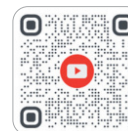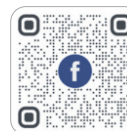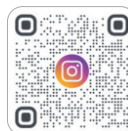

PRO ACTIV

PRO ACTIV Reha-Technik GmbH

Im Hofstätt 11
72359 Dotternhausen
info@proactiv-gmbh.de

www.proactiv-gmbh.com

Rev. 09/2026_01
Subject to technical modifications, printing errors & colour deviations.