



PROACTIV®



## WHEEL-E

### MAXIMUM DRIVE POWER AND CONTROL FOR ANY TERRAIN

The **WHEEL-E** towing device transforms your wheelchair into a powerful companion for any route. With its 20" wheel, wide tyres and low centre of gravity, you can master paved paths, forest trails and inclines with ease. The battery offers a range of up to 45 km, extendable to 90 km with a double battery. Thanks to cruise control and stable steering geometry, you can enjoy maximum comfort and safety. Versatile accessories enable customisation – for independent mobility on any terrain!




#### #IndependentOnTheRoad

Activate the throttle with a flick of the wrist and roll effortlessly over pavements, forest paths or inclines thanks to the 20" wheel and wide tyres.



#### #MoreSpeed

Depending on the equipment, you can reach speeds of 6 km/h, 15 km/h (optionally with individual operating licence) or 24 km/h. Whether in the city, off-road or on private roads – the choice is yours.

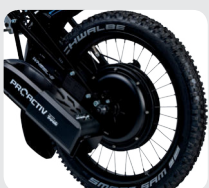

#### #SteepIncline

Up to 15 % climbing ability. With the powerful Triride hub drive, you can even tackle challenging inclines. Perfectly combined with 5 support levels and a reverse gear for more flexibility.



#### #DoubleSoFar

Whether you're travelling short or long distances: up to 45 km range with one battery. With the optional double battery, you can enjoy double the range – and the interchangeable plug connection makes it easy to change the battery.





## CUSTOMISED FOR YOU.

The **WHEEL-E** adapts perfectly to your needs: With an adjustable steering column in length, angle and height, you can find the ideal handlebar position. The flexible adapter can also be adjusted in height, length and angle – for maximum comfort and control.

## UNIQUE, SELECTABLE OPTIONS FROM PRO ACTIV:



The sturdy centre or front adapter can be fitted in a few simple steps and hardly adds any weight to the wheelchair.



Cruise control for relaxed journeys. Enjoy constant speed at the touch of a button. No adjustment, no stress – just pure driving pleasure.



Individuality counts. With accessories such as mudguards, pannier racks and rear-view mirrors, the WHEEL-E can be customised to your needs – even for functionally impaired users.



150 kg are already possible as standard. In individual cases, approval for higher load weights may be possible with special equipment. PRO ACTIV will check this on a case-by-case basis.






## TECHNICAL SPECIFICATIONS

Motor:	Hub motor: Triride I max. 1000 W
Battery:	Triride 48 V   10.5 Ah   504 Wh   Range approx. 25 - 45 km   3.3 kg
Double battery:	Possible; switch via change-over device on the bottom bracket support
Total weight without adapter:	Starting at 14.7 kg (without battery)
Max. standard load weight:	Up to 150 kg trailer load and 10 kg payload
Pack size:	Starting at 90-110 x 63 x 62 cm
Wheel sizes:	20"
Drive system speed:	6 km/h (Switzerland: 10 km/h), 15 km/h (only permitted with individual operating licence in Germany in public areas and in road traffic), 24 km/h (only permitted on private property in Germany)

Data subject to design changes.

## FOLLOW US ON SOCIAL MEDIA

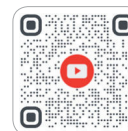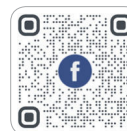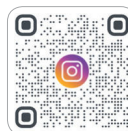

**PRO ACTIV®**

PRO ACTIV Reha-Technik GmbH

Im Hofstätt 11  
72359 Dotternhausen  
info@proactiv-gmbh.de

[www.proactiv-gmbh.com](http://www.proactiv-gmbh.com)

Rev. 05/2025  
Subject to technical modifications, printing errors & colour deviations.